

# READY HEAT

## 6 PANEL HEATED BLANKET

UP TO 10 HOURS OF HEAT



### INSTRUCTIONS FOR USE

1. DO NOT open the outer protective package until ready to use.
2. Remove the blanket from the outer protective package and expose the Ready-Heat blanket to air. DO NOT open blanket pockets or inside pouch.
3. Blanket will start to heat immediately. Allow 10-15 minutes to reach operating temperature.
4. Place the blanket on the patient with the Blue side Down, towards the patient and the Beige side Up, away from the patient.
5. If heat decreases, expose blanket to air.
6. After use, dispose of blanket with regular waste. All natural ingredients; will not harm environment.



### WARNINGS

\*\*Failure to follow these precautions may result in serious burns or other serious injuries.\*\*

1. **DO NOT use blanket on exposed skin.**
2. Careful supervision is required for use with the elderly, infants, children, people with sensitive skin and those who can't feel heat or regulate heat. Keep out of reach of children.
3. People with allergies, diabetes, frostbite, scars, open wounds or circulatory problems should consult a physician before using the Ready-Heat blanket.
4. Remove immediately if blanket becomes too hot or if you feel uncomfortable. Maximum temperature may reach 125°F (52°C). DO NOT use while sleeping.
5. DO NOT open the sealed pockets or pouch and DO NOT allow the contents to contact mouth or eyes. If contact with skin or eyes, wash with clean water. Contents contain iron powder which, if swallowed, can be harmful. If contents are swallowed, immediately call a poison control center.
6. DO NOT fold the blanket or place objects on top of heating elements as folding or pressure may increase heat.
7. DO NOT use in oxygen rich environment, including hyperbaric chamber or oxygen tent.

**CONTENTS:** Iron powder, water, activated carbon, wood powder and salt




## TechTrade




TechTrade, LLC ▪ 6900 Tavistock Lakes Blvd. ▪ Suite 400 ▪ Orlando, FL 32827 ▪ USA  
1-866-491-TECH ▪ Fax: 201-706-8131 ▪ info@techtradellc.com ▪ www.ready-heat.com

SIZE: 34" x 60" (86 x 152 cm)  
PART# S6RHLG  
WEIGHT: 2.5 lbs. (1.13 kg)  
NSN # 6532-01-525-4063

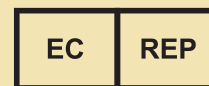

EMERGO EUROPE  
Prinsessegracht 20  
2514 AP The Hague  
The Netherlands

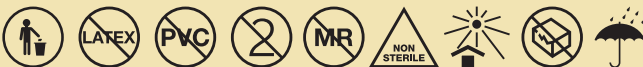

Please visit [www.ready-heat.com/IFU](http://www.ready-heat.com/IFU)  
for complete IFU and Symbol Descriptions.

MADE IN USA

0459

CLEARED



(01)00850017905023

RHPLG-2 rev 6, June 2021